BORMANN®

BWS1500

Art Nr: 025696

GB

GR

WWW.BORMANNTOOLS.COM

CE

Introduction

This product can be used for weighing applications as all normal weighing scales.

The device features a fat measuring system, and it can calculate weight, and at the same time, plan moisture rate measured fat rate, muscle a bone mass and basic metabolism etc. body composition.

Technical Data:

Available Modes: body fat, water, muscle, bone, bmi, calorie, weight

Max capacity:5kg-180kg/0.1kg

Measure fat range:5-50%/0.1%

Measure moisture range:10-8

Memory data: 12 profiles

high precision sensors,auto on/off

Overload Indication: err

Low Power indication: lo

Available Units: kg, lb, st

Graduation: 100g/0.2lb

Max. Capacity: 180kg,396lb

Min capacity: 5kg

Tempered safety glass platform: 280*280*5mm

Display: Lcd or led

Size of display:75*30mm

Battery Type:2*AAA

Operating Instructions

Selecting Personal Parameters

1. Press **SET** Key to turn the device on (**every time turn on scale to set it**). Select one of the flashing personal data results/ profiles (n=0-11), press the up arrow key or the down arrow key (to increase or decrease and choose the number which you want (the device can store up to 12 settings) then press **SET** to confirm. If there is no action within 3 seconds, system will confirm automatically the current flashing number is the current user measurement data.
2. After storing the user profile, the gender selection symbol will start flashing. Press the up arrow key or the down arrow key to select your gender. Press **SET** to confirm, if there is no operation, it will turn off automatically within 20 seconds.
3. After storing the gender, the age symbol will start flashing. Press the up arrow key or the down arrow key to adjust the age value, press **SET**.
To conform if there is no operation, it will turn off automatically within 20 seconds.
4. After storing your age, the height symbol will start flashing. Press the up arrow key or the down arrow key to adjust the height value, press **SET**.
To conform if there is no operation, it will turn off automatically within 20 seconds.
5. After storing all the data, the device will display "0.0", which indicates that you can proceed to the calculation of fat.

Measurements

A. Calculating weight

- 1) Under the turn off, please move on repose fully the scale, screen will show your weight, after 3 times, it will be locked.
- 2) If there is no operation,it will be turn off automatically within 10 seconds.

B. Measurement: Body fat, water, muscle, bone, calorie, weight, BMI

Please bare feet stand steadily on the scale ,feet with two electrode contact, body don't swing. The screen will show your weight value, stability weight data flashing three time then be locked, it will indicate fat,water, muscle, bone, calorie, weight ,each measurement value will be displayed for two seconds, alternating cycle three times after the automatic shut-down.


Usage Notes

1. Ensure that the scale is placed on a hard & level surface, otherwise the result will be inaccurate if on the carpet or uneven ground.
2. Keep scale surface dry to avoid accident by slipping.
3. When weighing please take off the socks, shoes and clothes, at relatively fixed time before meal.
4. In order to measure precision , not recommend to weight in 2 hours after getting up.
5. Please put the scale in a dry place to protect the internal electronic components.
6. Please use micro cloth to clean the scale, not immersed in water.
7. Please don't remove the scale except replace the battery.
8. Please don't use heavy percussion or falls from a height.
9. Please don't persistent put heavy on the scale to avoid damage to the high precision sensor.

* The manufacturer reserves the right to make minor changes to product design and technical specifications without prior notice unless these changes significantly affect the performance and safety of the products. The parts described / illustrated in the pages of the manual that you hold in your hands may also concern other models of the manufacturer's product line with similar features and may not be included in the product you just acquired.

* Please note that our equipment has not been designed for use in commercial, trade or industrial applications. Our warranty will be voided if the machine is used in commercial, trade or industrial businesses or for equivalent purposes.

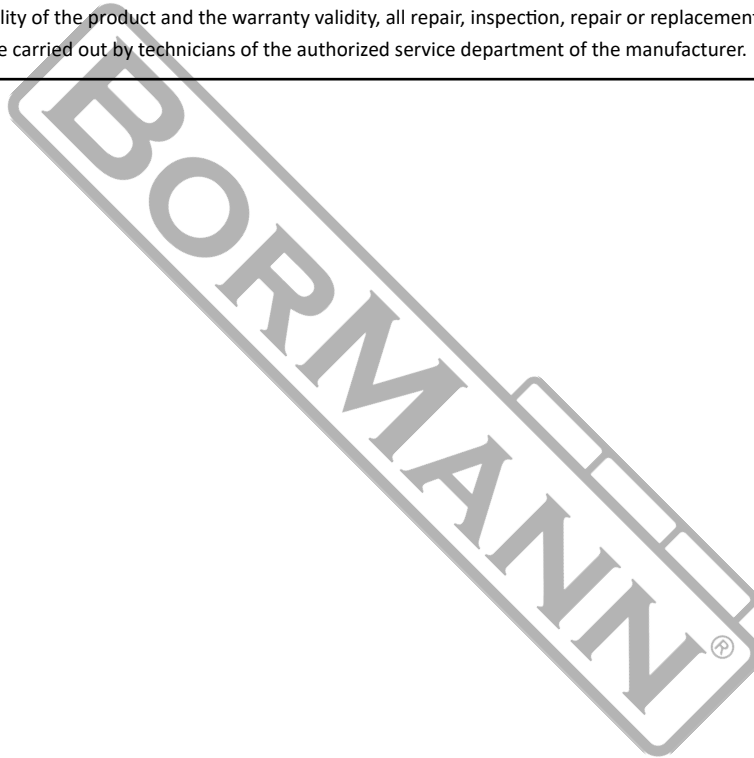* To ensure the safety and reliability of the product and the warranty validity, all repair, inspection, repair or replacement work, including maintenance and special adjustments, must only be carried out by technicians of the authorized service department of the manufacturer.

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